

Neo Cot

Bassinet



Neo Cot

A bassinet to bring mother and baby closer together.

A bassinet to encourage interaction between mother and baby, which promotes eye contact and visibility of the baby.



The transparent cot can be locked.

The transparent cot can be locked at 4 points underneath the mattress.



“Interaction”

The height of the HL stand can be adjusted to promote eye contact between mother and baby.

The lowest position of the HL stand is 62cm which enables the cot to be adjusted to bedside height.

The Neo Cot has been designed to support the caring model of both mothers and nursing staff, paying special attention to safety and interaction.

Caster-lock from any direction

Four-wheel locking mechanism ensures the cot remains stable when interacting with the baby.



Mattress can be tilted using a simple lever situated at the front of the cot

Stepless smooth tilting is possible while watching a baby.



Easy to move

Large dual wheel casters make it easy to operate and move the cot.



A combination of options can be selected to meet a variety of nursery requirements.

Enhances operation efficiency for nursing staff

Variety of peripheral devices can be attached to the IV Pole (option), which enables a staff to operate the devices immediately in a large work space.

Drawer

Pulse oximeter can be put on the slide tray which is suitable for practical use.



Sliding open the top slide tray


Opening the drawer

IV Pole

IV Pole can be attached to each corner.

Basket

Small items such as disinfectant and thermometers can be stored.

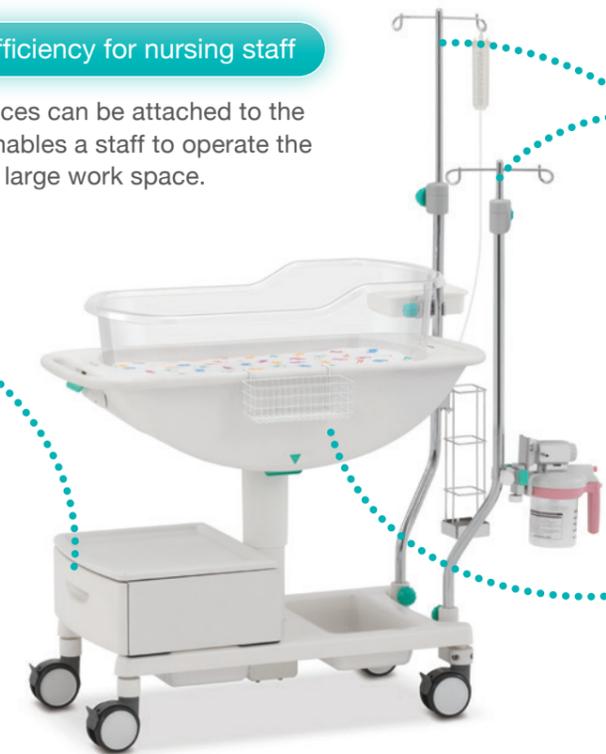

Carefully selected material ensures improved ease of cleaning.

Increased ease of cleaning of the product

The main body of the Neo Cot is made of resin which allows for a variety of different cleaning materials to be used.

*Please confirm supplier of cleaning material if it could be used for resin.



The mattress is washable

Highly breathable and washable material.



Specification

		Neo Cot HL Stand
Atom's Code		26408
Tilting Angle		0° - 12°
Mattress	Material	Polyester-elastomer
	Method of disinfection	Dipping (avoid using dry heat sterilization)
Caster		Dia. 100 mm (Dual caster w/stopper)
Load Capacity		10 kg (on mattress)
Weight		Approx. 27 kg
Dimensions		W960 x D450 x H620 - 820 mm

Options & Supplies



Atom's code: 26113
I.V Pole for Neo Cot



Atom's code: 26114
Drawer for Neo Cot



Atom's code: 26117
Basket for Neo Cot



Atom's code: 26411
Cylinder Rack for Neo Cot
 Description: Users can use cylinders with the outer diameter of ≤105 mm.



Atom's code: 26123
Shade for Neo Cot




Atom's code: 26118
Mattress for Neo Cot



Atom's code: 26122
Mattress Sheet for Neo Cot (5 pcs/pkg)

- Atom's code: 26112 **Transparent Cot for Neo Cot**
- Atom's code: 60337 **Neo Sheets (100 pcs/box)**
- Atom's code: 67100 **Baby card , blue (50 pcs/pkg)**
- Atom's code: 67101 **Baby card , pink (50 pcs/pkg)**
- Atom's code: 26124 **Basket wide type for Neo Cot**

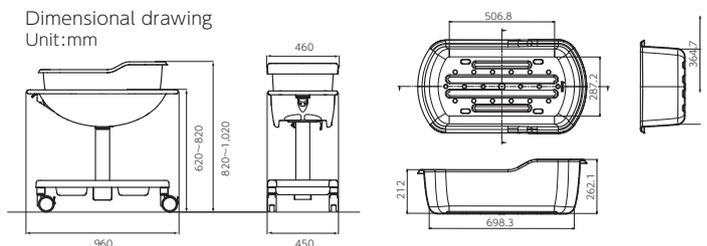

Atom Medical products are produced under strict production control in Japan.
 We provide value and trust of the highest quality for the tiny baby's lives in the world.